HARTMANN

Atrauman® Ag

Silberhaltige Salbenkompress
Silver-containing impregnated dressing
Comresse d'argent imprégnée d'un excipient
Zilverhoudend zalfkompres
Compressa con pomata contenente argento
Apósito de plata impregnado con pomada
neutra
Compressa não aderente com prata
Επίθεμα με αλοιφή και προσθήκη αργύρου
Kompres s mastí obsahující stříbro
Kompres s mast'ou obsahující striebro
Opatrunek z maścią zawierający srebro
Ezüst tartalmú, kenőccsel impregnált kötés
Серебросодержащая мазевая повязка
Silverhaltig salvkompres
Sølvholdigt salvekompres
Sølvholdig salvekompres
含銀軟膏敷料
ضمادة ذات مرهم يحتوي على الفضة

STERILE

R

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Popis výrobku

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Instructions for use

Product description

Atrauman Ag is a non adherent tulle dressing containing silver with anti-bacterial properties for atraumatic wound treatment.

Composition

Support fabric made of polyamide fibres coated with elemental silver. The impregnation is non-medicated and triglyceride based (neutral lipids).

The impregnation contains (INCI assigned names): Caprylic/Capric/Stearic Triglyceride; Bis-diglyceryl-polyacyladipate-2 and Macrogol 2000.

Properties and modes of action

The soft, thin support fabric drapes easily and ensures close contact with the wound base. The fabric's metallic surface and the impregnation counteract adhesion to the wound; as a result the dressing can be removed painlessly.

Bacteria on the dressing surface are destroyed by the silver ions created in the moist environment. Atrauman Ag works against Gram-negative (e.g. *Klebsiella pneumoniae*) as well as Gram-positive bacteria (e.g. *Staphylococcus aureus* including the MRSA strains).

The impregnation protects the wound edges helping to prevent maceration. It is non-medicated, contains no paraffin and leaves no residues. Because of this, when the dressing is changed there is no residue left in the wound or on the wound edges.

Indications

Atrauman Ag is suitable for the treatment of superficial, acute and chronic wounds of all kinds; it is particularly suitable for complementary use in the treatment of infected wounds or wounds that show a high bacterial contamination. Atrauman Ag is also indicated as an additional infection prophylaxis in wound treatment.

Mode of application

Take the Atrauman Ag from the peel pack with both cover papers in place, remove one of the papers and place this side on the wound. Remove the remaining cover paper.

Secure a sterile, absorbent dressing pad e.g. Zetuvit, over Atrauman Ag to absorb the exudate. The microbiological action of Atrauman Ag lasts for at least seven days. Unless otherwise prescribed by the doctor or medically indicated, a new Atrauman Ag dressing can be applied at each dressing change.

In case of thick exudate it is recommended to cut Atrauman Ag to prevent congestion of secretions.

Contraindications

None known.

Specific precautions

Before treating any infected wounds a medical assessment of the wound condition and the causes of infection must first be obtained. Treatment with Atrauman Ag cannot replace the requisite antibiotic therapy.

So as to not reduce the efficacy of Atrauman Ag, the wound dressing should not be used in combination with dressings or ointments containing paraffin.

Atrauman Ag must not be used in combination with iodine-containing preparations.

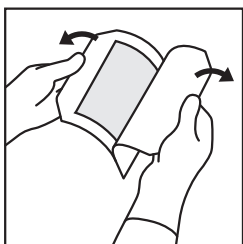
Special instructions

Store dry in a horizontal position at room temperature.

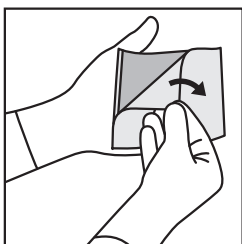
Contents sterile unless peel packs are damaged.

Do not resterilise. Keep out of the reach of children.

Date of revision of the text: 2009-04

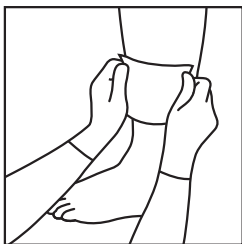
1

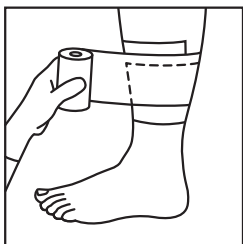
2

3

4

5